



*Ahsanullah University of Science and Technology*

*Department of Arts and Sciences*

**LABORATORY MANUAL**  
**FOR**  
**CHEMISTRY SESSIONAL COURSES**

**Manual I**  
**Quantitative Inorganic Analysis**

# Chemistry Manual I

## Inorganic Quantitative Analysis

AUST

### Experiment No. 01

**Experiment Name:** Standardization of sodium hydroxide solution with a standard solution of oxalic acid.


#### Preparation of standard 0.050 M oxalic acid solution

Weigh out exactly 0.63 g of pure crystalline oxalic acid dihydrate ( $\text{HOOC-COOH}\cdot 2\text{H}_2\text{O}$ ) in a 100 ml volumetric flask. Dissolve it with little amount of distilled water and fill it up to the mark.

#### Working Procedure

- i) Take the supplied sodium hydroxide solution in the burette.
- ii) Pipette out exactly 10.0 ml of the prepared oxalic acid solution into a conical flask.
- iii) Add 1-2 drops of phenolphthalein as indicator.
- iv) Record the initial burette reading.
- v) Add sodium hydroxide solution from the burette into the conical flask with constant shaking until the color changes from colorless to pink and record the burette (final reading) again.
- vi) Repeat the procedures (ii-v) at least three times and take the average (for average, data which have differences within 0.3 should be considered). This value is the volume of sodium hydroxide solution required to neutralize the standard 10.0 ml oxalic acid solution.
- vii) Calculate the concentration in terms of Molarity by using both formula and unitary method.

#### Reaction

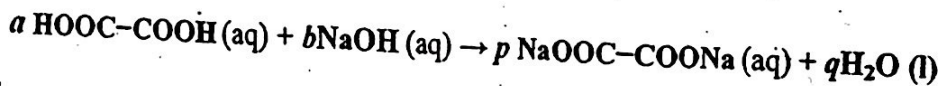

#### Calculation:

From the above reaction, it is clear that in order to be neutralized;  $a$  mol (coefficient) of oxalic acid reacts with  $b$  mol (coefficient) of sodium hydroxide. The concentration of sodium hydroxide solution (in terms of molarity) can be determined using the following formula:

$$\frac{M_1 V_1}{a} = \frac{M_2 V_2}{b}$$

Where  $M_1$  = Molarity of oxalic acid solution = 0.050 M

$V_1$  = Volume of oxalic acid solution taken = 10.0 ml

$M_2$  = Molarity of sodium hydroxide solution = ? M

$V_2$  = Volume of sodium hydroxide solution required = Avg. burette reading (? ml)

#### Questions:

1. Which one of NaOH and Oxalic acid is the primary standard substance in this experiment? Explain giving definition of primary standard and secondary standard substance.
2. Why Phenolphthalein is chosen as the suitable indicator for this experiment? Use titration curve for explanation.
3. Write the balanced reaction. Mention the type of reaction with definition.
4. Calculate the strength of NaOH solution in normality (N).

Fi/

Experiment No. 01

**Experiment Name:** Standardization of sodium hydroxide solution with a standard solution of oxalic acid.


**Preparation of standard 0.050 M oxalic acid solution**

Weigh out exactly **0.63 g** of pure crystalline oxalic acid dihydrate (HOOC-COOH.2H<sub>2</sub>O) in a **100 ml** volumetric flask. Dissolve it with little amount of distilled water and fill it up to the mark.

**Working Procedure**

- i) Take the supplied sodium hydroxide solution in the burette.
- ii) Pipette out exactly **10.0 ml** of the prepared oxalic acid solution into a conical flask.
- iii) Add **1-2 drops of phenolphthalein** as indicator.
- iv) Record the initial burette reading.
- v) Add **sodium hydroxide** solution from the burette into the conical flask with constant shaking until the color changes from **colorless** to **pink** and record the burette (final reading) again.
- vi) Repeat the procedures (ii-v) at least three times and take the average (for *average, data which have differences within 0.3 should be considered*). This value is the volume of sodium hydroxide solution required to neutralize the standard 10.0 ml oxalic acid solution.
- vii) Calculate the concentration in terms of **Molarity** by using both **formula** and **unitary method**.

**Reaction**

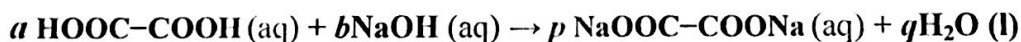

**Calculation:**

From the above reaction, it is clear that in order to be neutralized; **a mol (coefficient)** of oxalic acid reacts with **b mol (coefficient)** of sodium hydroxide. The concentration of sodium hydroxide solution (in terms of molarity) can be determined using the following formula:

$$\frac{M_1 V_1}{a} = \frac{M_2 V_2}{b}$$

Where  $M_1$  = Molarity of oxalic acid solution = 0.050 M

$V_1$  = Volume of oxalic acid solution taken = 10.0 ml

$M_2$  = Molarity of sodium hydroxide solution = ? M

$V_2$  = Volume of sodium hydroxide solution required = Avg. burette reading (? ml)

**Questions:**

1. Which one of NaOH and Oxalic acid is the primary standard substance in this experiment? Explain giving definition of primary standard and secondary standard substance.
2. Why Phenolphthalein is chosen as the suitable indicator for this experiment? Use titration curve for explanation.
3. Write the balanced reaction. Mention the type of reaction with definition.
4. Calculate the strength of NaOH solution in normality (N).



# Chemistry Manual I

## Inorganic Quantitative Analysis

AUST

### Experiment No. 02


**Experiment Name:** Standardization of HCl solution with a secondary standard solution of sodium hydroxide (NaOH).

**Standardization of sodium hydroxide solution:** The supplied sodium hydroxide solution was standardized as per the procedure described in **exp # 01**

#### Working Procedure

- i) Take the standard sodium hydroxide (NaOH) solution in the burette. <sup>(a)</sup>
- ii) Pipette out exactly **10.0 ml** of the supplied hydrochloric acid solution into a conical flask by a volumetric pipette.
- iii) Add **1-2 drops of methyl orange** as an indicator.
- iv) Record the initial burette reading.
- v) Titrate the hydrochloric acid solution with the sodium hydroxide solution until the color changes from **pink to yellow**. <sup>(b)</sup>
- vi) Note the final burette reading.
- vii) Repeat the procedure (ii-vi) at least three times and take the average. This value is the volume of sodium hydroxide solution required to neutralize the hydrochloric acid solution.
- viii) Calculate the concentration in terms of **Molarity** by using both **formula and unitary methods**.

#### Reaction:

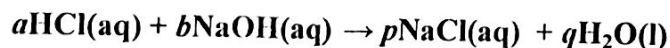

#### Calculation:

Calculate the molarity of the HCl solution from the following formula:

$$M_1V_1 / a = M_2V_2 / b$$

Where  $M_2$  = Molarity of sodium hydroxide solution.

$V_2$  = Volume of sodium hydroxide solution taken.

$M_1$  = Molarity of hydrochloric acid solution (=?).

$V_1$  = Volume of hydrochloric acid solution required.

#### Questions:

1. Show by pH – Neutralization curve that methyl orange is suitable for this titration.
2. Do you think that phenolphthalein or methyl red can be used here? Justify your answer.
3. Calculate the concentration of HCl in terms of normality (N).

#### Note:

- (a) If this experiment is not carried out just after finishing the experiment-01 then HCl solution is better to take in the burette instead of NaOH solution.
- (b) In some places color change is written as red to orange but pink to yellow is usually observed.



Experiment No. 02


Experiment Name: Standardization of HCl solution with a secondary standard solution of sodium hydroxide (NaOH).

Standardization of sodium hydroxide solution: The supplied sodium hydroxide solution was standardized as per the procedure described in exp # 01

Working Procedure

- i) Take the standard sodium hydroxide (NaOH) solution in the burette. <sup>(a)</sup>
- ii) Pipette out exactly 10.0 ml of the supplied hydrochloric acid solution into a conical flask by a volumetric pipette.
- iii) Add 1-2 drops of methyl orange as an indicator.
- iv) Record the initial burette reading.
- v) Titrate the hydrochloric acid solution with the sodium hydroxide solution until the color changes from pink to yellow. <sup>(b)</sup>
- vi) Note the final burette reading.
- vii) Repeat the procedure (ii-vi) at least three times and take the average. This value is the volume of sodium hydroxide solution required to neutralize the hydrochloric acid solution.
- viii) Calculate the concentration in terms of Molarity by using both formula and unitary methods.

Reaction:

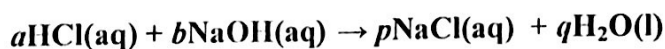

Calculation:

Calculate the molarity of the HCl solution from the following formula:

$$M_1V_1 / a = M_2V_2 / b$$

Where  $M_2$  = Molarity of sodium hydroxide solution.

$V_2$  = Volume of sodium hydroxide solution taken.

$M_1$  = Molarity of hydrochloric acid solution (=?).

$V_1$  = Volume of hydrochloric acid solution required.

Questions:

1. Show by pH – Neutralization curve that methyl orange is suitable for this titration.
2. Do you think that phenolphthalein or methyl red can be used here? Justify your answer.
3. Calculate the concentration of HCl in terms of normality (N).

Note:

- (a) If this experiment is not carried out just after finishing the experiment-01 then HCl solution is better to take in the burette instead of NaOH solution.
- (b) In some places color change is written as red to orange but pink to yellow is usually observed.

Experiment No. 03

**Experiment Name:** Standardization of hydrochloric acid (HCl) solution with a standard solution of sodium carbonate (Na<sub>2</sub>CO<sub>3</sub>).


**Preparation of standard Na<sub>2</sub>CO<sub>3</sub> solution:**

Weigh out exactly 0.53 g of pure crystalline sodium carbonate (Na<sub>2</sub>CO<sub>3</sub>) in a 100 ml volumetric flask. Dissolve it with a little amount of distilled water and fill it up to the mark.

**Working Procedure:**

- i) Take the supplied hydrochloric acid (precaution must be maintained) solution in the burette.
- ii) Pipette out exactly 10.0 ml of the prepared sodium carbonate solution into a conical flask by a volumetric pipette.
- iii) Add 1-2 drops of methyl orange as an indicator.
- iv) Record the initial burette reading.
- v) Titrate the sodium carbonate solution with the hydrochloric acid solution until the first drop changes its color from yellow to pink.
- vi) Note the final burette reading.
- vii) Repeat the procedure (ii-vi) at least three times and take the average. This value is the volume of hydrochloric acid solution required to neutralize the Na<sub>2</sub>CO<sub>3</sub> solution.
- viii) Calculate the concentration in terms of Molarity by using both formula and unitary methods.

**Reaction:**

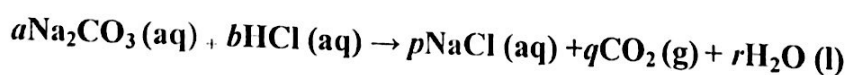

**Calculation:** Calculate the molarity of the hydrochloric acid solution from the following formula:

$$M_1V_1 / a = M_2V_2 / b$$

Where, M<sub>1</sub> = Molarity of sodium carbonate solution.  
V<sub>1</sub> = Volume of sodium carbonate solution taken.  
M<sub>2</sub> = Molarity of hydrochloric acid solution.  
V<sub>2</sub> = Volume of hydrochloric acid solution.

**Questions:**

1. The above reaction actually occurs in two steps. Show these steps.
2. What is salt? What type of salt sodium carbonate is: acidic or basic? Explain with reaction.
3. Sketch the titration curve and give the reason for selecting methyl orange as indicator.
4. Can you explain why indicator changes its color? **Hint:** theories of indicator action.
5. Calculate the strength of HCl in normality (N).



Experiment No. 04

**Experiment Name:** Standardization of sodium thiosulphate ( $\text{Na}_2\text{S}_2\text{O}_3$ ) solution with a standard solution of potassium dichromate ( $\text{K}_2\text{Cr}_2\text{O}_7$ ).



**Preparation of standard potassium dichromate ( $\text{K}_2\text{Cr}_2\text{O}_7$ ) solution**

Weigh out exactly **0.47 g** of potassium dichromate ( $\text{K}_2\text{Cr}_2\text{O}_7$ ) in a **100 mL** volumetric flask. Dissolve it with a little volume of distilled water by shaking. Make the solution up to the mark by adding distilled water.

**Working Procedure**

- i) Take the supplied sodium thiosulphate ( $\text{Na}_2\text{S}_2\text{O}_3$ ) solution in the burette & note the initial reading.
- ii) Take about **50 ml** of distilled water into a conical flask by using measuring cylinder.
- iii) Add about **1-2 g** of potassium iodide (**KI**) into this conical flask.
- iv) Add about **1 g** of sodium bicarbonate and shake well to dissolve it.
- v) Add about **3 ml** of **conc. hydrochloric acid** slowly while gently rotating the flask to mix the liquids. ( *$\text{CO}_2$  evolved here to drive out  $\text{O}_2$  from solution and to blanket/cover the solution surface. Check the solution by a piece of litmus paper whether it is acidic or not.*)
- vi) Add exactly **10.0 ml** of potassium dichromate solution and mix the solution well.
- vii) Cover the flask immediately with a watch glass and allow standing in a **dark box** (*to prevent the solution from light*) for **5 minutes**. (*solution will be dark brown*).
- viii) Take out the flask from the dark box and dilute the solution by adding **50 mL of distilled water**.
- ix) Titrate the liberated iodine with the sodium thiosulphate solution from the burette until the color **fades and turns pale yellow**. (compare with  $\text{K}_2\text{Cr}_2\text{O}_7$  solution)
- x) Now add **1 ml or 10 drops of starch** (*indicator; indicates iodine*) solution into it. The solution will become **deep violet**.
- xi) Continue the addition of sodium thiosulphate solution until the solution changes its color to **light green**.
- xii) Note the final burette reading. This is the volume of sodium thiosulphate solution required.
- xiii) Repeat the procedure (ii-xii) at least three times and take the average value.
- xiv) Calculate the concentration in terms of **Molarity** by using both **formula** and **unitary methods**.

**Reactions:**

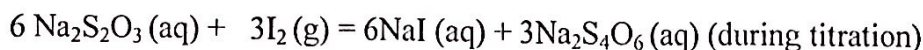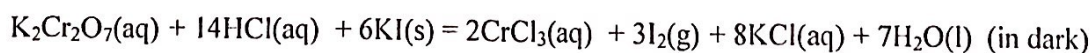

Therefore,  $1 \text{ mol K}_2\text{Cr}_2\text{O}_7 \equiv 3 \text{ mol I}_2 \equiv 6 \text{ mol Na}_2\text{S}_2\text{O}_3$

And finally,  $1 \text{ mol K}_2\text{Cr}_2\text{O}_7 \equiv 6 \text{ mol Na}_2\text{S}_2\text{O}_3$



**Calculation:**

Calculate the molarity of the sodium thiosulphate solution from the following formula:

$$M_1V_1 / 1 = M_2V_2 / 6$$

Where  $M_1$  = Molarity of potassium dichromate solution

$V_1$  = Volume of potassium dichromate solution taken

$M_2$  = Molarity of sodium thiosulphate solution (= ? M)

$V_2$  = Volume of sodium thiosulphate solution used.

**Questions:**

1. Define oxidation and reduction reactions in terms of electron transferred and in terms of oxidation number?
2. Give half-reactions for  $K_2Cr_2O_7$ , KI,  $Na_2S_2O_3$  and  $I_2$  and then level them as oxidizing and reducing agents with explanation?
3. What are Iodometric and Iodimetric titrations? Give example.
4. What is the function of  $NaHCO_3$  used in this experiment? Why is the solution kept in the dark box/place and for at least 5 minutes?
5. Write the structure of starch. Mention the advantages of using starch as an indicator in this experiment. Indicator is used in this experiment at the last stage – why?

Experiment No. 05

**Experiment Name:** Estimation of the amount of copper present in a supplied solution by iodometric method.




**Standardization of the sodium thiosulphate solution**

The supplied sodium thiosulphate solution was standardized as per procedure described in (Exp. 4).

**Working Procedure:**


- (i) Take the supplied sodium thiosulphate solution into the burette.
- (ii) Take exactly 10.0 ml of the supplied copper solution with a volumetric pipette into a conical flask.
- (iii) Add one drop of concentrated acetic acid solution into the conical flask.
- (iv) Add about 1 g of solid potassium iodide (KI) into it and cover it with a watch glass.
- (v) Keep the flask in a dark box for 5 minutes. A reddish brown color will be formed. Then add 20 ml distilled water into conical flask.
- (vi) Titrate the solution with standard sodium thiosulphate solution until the color becomes light yellow. (*this is not the end point*)
- (vii) Add 10 drops of starch solution into the conical flask.
- (viii) Add about 10 ml of 10% NH<sub>4</sub>CNS solution and continue adding sodium thiosulphate until the color of the solution becomes straw/pale (white). (*it is the end point*)
- (ix) Note the final burette reading.
- (x) Repeat the titration at least three times and take the average value.

**Reactions:**

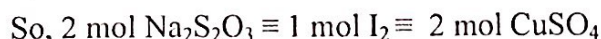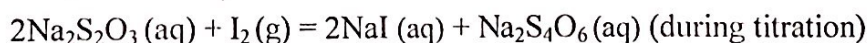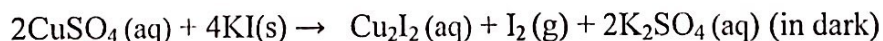

**Calculation:**

**Method 1:** Calculate the molarity of CuSO<sub>4</sub> solution. Now use the following relationship:

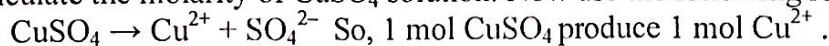

Therefore, the concentration of CuSO<sub>4</sub> solution = the concentration of Cu<sup>2+</sup>

$$\text{Now, } \frac{1 \text{ M Cu}^{2+}}{\text{mole}} = 1 \text{ mol/L Cu}^{2+} = 63.54 \text{ g/L Cu}^{2+}$$

**Method 2:** 1M thiosulphate solution  $\equiv$  1M CuSO<sub>4</sub> = 1 mol/L CuSO<sub>4</sub>  
= 1 mol/L Cu<sup>2+</sup> = 63.54 g/L Cu<sup>2+</sup>

So, 1000 ml of 1M thiosulphate solution = 63.54 g of Cu<sup>2+</sup>

And finally, 1ml of 1M thiosulphate solution = 0.06354 g Copper.

**Questions:**

1. Show the half-reactions involved here.
2. What is the function of NH<sub>4</sub>CNS in this titration?



# Chemistry Manual I

## Inorganic Quantitative Analysis

AUST

### Experiment No. 06

**Experiment Name:** Standardization of potassium permanganate solution ( $\text{KMnO}_4$ ) with a standard sodium oxalate ( $\text{Na}_2\text{C}_2\text{O}_4$ ) solution.

#### Preparation of Standard sodium oxalate solution:

Weigh out accurately **0.67 g** of sodium oxalate salt in a **100 ml** volumetric flask. Dissolve by adding a little amount of distilled water and make up to the mark with distilled water.

#### Working Procedure

- (i) Take the supplied potassium permanganate ( $\text{KMnO}_4$ ) solution into the burette and note the initial reading.
- (ii) Pipette out exactly **10.0 ml** of sodium oxalate ( $\text{Na}_2\text{C}_2\text{O}_4$ ) solution by a volumetric pipette into a conical flask.
- (iii) Add about **50 ml** of **1M Sulphuric acid** solution to the conical flask.
- (iv) Carry out the titration (add 2 – 5 drops of  $\text{KMnO}_4$  solution) at the ordinary temperature until the **first point pink** color appears throughout the solution. Allow to stand until the pink color disappears and the solution becomes **colorless**.
- (v) Warm the solution to **50-60 °C** and continue titration to a **first permanent pink color**.
- (vi) Repeat the titration at least three times and take the average value.
- (vii) Calculate the concentration in terms of **Molarity** by using both **formula** and **unitary methods**.

#### Reaction:



#### Calculation

$$M_1 V_1 / \cancel{b} = M_2 V_2 / \cancel{a}$$

where  $M_1$  = Molarity of sodium oxalate solution  
 $V_1$  = Volume of sodium oxalate solution  
 $M_2$  = Molarity of potassium permanganate solution (=?)  
 $V_2$  = Volume of potassium permanganate solution

#### Questions:

1. Balance the above equation by half-reaction method.
2. Label  $\text{Na}_2\text{C}_2\text{O}_4$  and  $\text{KMnO}_4$  as oxidizing and reducing agents with reaction.
3. Sulphuric acid is used here; hydrochloric acid or nitric acid is not used. Explain why?
4. Find the normality of  $\text{Na}_2\text{C}_2\text{O}_4$  and  $\text{KMnO}_4$  solutions in this experiment.
5. No indicator was used in this experiment. Why?



Experiment No. 07

**Experiment Name:** Estimation of ferrous ions ( $\text{Fe}^{2+}$ ) present in a supplied solution of ferrous sulphate ( $\text{FeSO}_4$ ) by standard potassium permanganate ( $\text{KMnO}_4$ ) solution.


Standardization of the supplied potassium permanganate solution

Standardize the potassium permanganate solution by the following the procedure described in (expt.6).

Working Procedure

- (i) Take the standard potassium permanganate solution into the burette.
- (ii) Pipette out **exactly 10.0 ml** of supplied ferrous sulphate solution by a volumetric pipette into the conical flask.
- (iii) Add about **10 ml of 1M sulphuric acid** solution.
- (iv) Titrate the solution with a standard potassium permanganate ( $\text{KMnO}_4$ ) solution to the first **permanent pink color**.
- (v) Repeat the titration at least three times and take the average value.

Reaction:

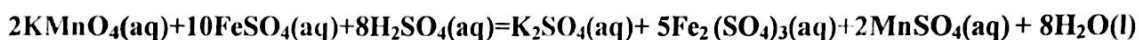

Therefore,  $2 \text{ mol KMnO}_4 \equiv 10 \text{ mol FeSO}_4$

or,  $1 \text{ mol KMnO}_4 = 5 \text{ mol FeSO}_4$



Calculation:

Method-1:

Calculate the molarity of  $\text{FeSO}_4$  solution. To find the amount of  $\text{Fe}^{2+}$  per litre solution, see the calculation of **Method 1** in Experiment #5.

Here,  $1 \text{ mol Fe}^{2+} = 55.85 \text{ g Fe}^{2+}$

Method-2: Calculate the amount of iron per liter of solution from the following equations:

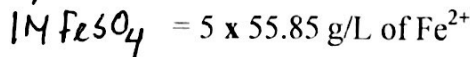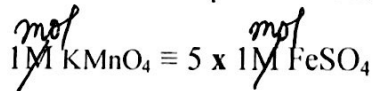

So, ~~1 ml of 1M KMnO<sub>4</sub> = 5 x 0.05585 g of Fe<sup>2+</sup>~~ (see Method 2 in Exp. 5)

Questions:

1. Make the above equation by using half-reactions.
2. What is auto indicator? Give two examples.

\* Auto catalyst, \* Self indicator, \* Temp. effect on reaction kinetics.

Experiment No. 08

**Experiment Name:** Determination of total Hardness of water using Eriochrome Black T (EBT) as an indicator.

**Hardness**

Water hardness was understood to be a measure of the capacity of water to precipitate soap. Soap is precipitated chiefly by calcium and magnesium ions present. Other polyvalent cations also may precipitate soap but they often are in complex forms, frequently with organic constituents, and their role in water hardness may be minimal and difficult to define. Total hardness is defined as the sum of the calcium and magnesium concentration, both expressed as calcium carbonate in mg/L (milligram/litre). The hardness may range from zero to hundreds of milligrams per litre, depending on the source and treatment to which has been subjected.

**Basic Principle**

Ethylenediaminetetraacetic acid (EDTA) forms a chelated soluble complex when added to a solution of certain metal cations. If a small amount of dye such as Eriochrome Black T (EBT) is added to an aqueous solution containing calcium and magnesium ions at a  $P^H$  of  $10.0 \pm 0.1$ , the solution becomes wine red. If EDTA is added as a titrant, the calcium and magnesium will be complexed, and when all the magnesium and calcium has been complexed the solution turns from wine red to blue, making the end point of the titration.

**Preparation of 0.01M EDTA disodium salt (EDTA) solution:**

Weigh out accurately **0.372 g of EDTA ( $C_{10}N_2O_8H_{14}Na_2 \cdot 2H_2O$ )** powder into a **100ml** volumetric flask and add a little water to dissolve it. Make up to mark with distilled water.

**Working Procedure**

- (i) Take the **standard EDTA** solution into the burette.
- (ii) Take **50 ml of sample water** (tap water) by a volumetric pipette into a conical flask.
- (iii) **Acidify** the sample water with **conc. hydrochloric acid** (use a litmus paper) and boil for a minute to drive off carbon dioxide ( $CO_2$ ).
- (iv) **Cool** and neutralize it with sodium hydroxide (use litmus paper).
- (v) Add about **2 ml of buffer solution** (pH 10) and **2 drops of EBT** indicator.
- (vi) Titrate the solution with a standard **0.01M EDTA** solution until the color changes from **wine red to blue**
- (vii) Repeat the titration at least three times and take the mean value. The mean of three readings of EDTA is the volume of EDTA required to calculate the total hardness.

**Calculation**

**Method 1** Unitary method

Find the concentration of  $\text{Ca}^{2+}$  using 1:1 complex formation with EDTA.  
Then use the relationship between  $\text{Ca}^{2+}$  and  $\text{CaCO}_3$  to calculate the amount (mg) of  $\text{CaCO}_3$  per litre of sample water. (*Be careful about the volume of water; it is not 10 ml*)

**Method 2** (see hints below)

1 ml 0.01 M EDTA solution = 1.00 mg of  $\text{CaCO}_3$  (show by calculation)

Therefore, total hardness as mg  $\text{CaCO}_3$  per liter of water

$$= (a \times b \times 1000) \div \text{ml of sample taken}$$

Where, a = ml. of titrant required (mean volume of EDTA solution)


b = mg  $\text{CaCO}_3$  equivalent to 1.00 ml EDTA titrant. (100 mg)

**Questions:**

1. Discuss on hard water and the process to soften it.
2. What are the permanent and temporary hardness of water?
3. Define the terms Chelete and polydentate.
4. Write down the structure of EDTA. What type of reaction takes place between the Calcium and EDTA?
5. What is soap? What is meant by micelle? How soap is related with hard water or hardness of water?

**Hints:**

[Molar mass of  $\text{CaCO}_3$  is 100 g/mol]

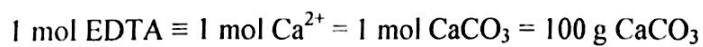

So 1M EDTA  $\equiv$  1M  $\text{CaCO}_3$  = 1mol/L  $\text{CaCO}_3$  = 100 g/L  $\text{CaCO}_3$

Hence 1000 ml 0.01M EDTA = 0.01 x 100 g  $\text{CaCO}_3$

$$= 0.01 \times 100 \times 1000 \text{ mg CaCO}_3$$

$$= 1000 \text{ mg CaCO}_3$$



Experiment No. 09

**Name of the Experiment: Determination of the available Chlorine in bleaching Powder.**


**Theory:**

Bleaching powder consists essentially of a mixture of calcium hypochlorite,  $\text{Ca}(\text{OCl})_2$  and the basic calcium chloride  $\text{CaCl}_2 \cdot \text{Ca}(\text{OH})_2 \cdot \text{H}_2\text{O}$ . Some free slaked lime is usually present. The active constituent is the hypochlorite ( $\text{OCl}^-$ ), which is responsible for bleaching action. Upon treating bleaching powder with hydrochloric acid, chlorine is liberated.



The available chlorine refers to the chlorine liberated by the action of dilute acid, and is expressed as the percentage by weight of the bleaching powder.

The available chlorine in bleaching powder can be determined iodometrically. For this the bleaching powder solution or suspension is treated with an excess of a solution of potassium iodide, and strongly acidified with acetic acid.

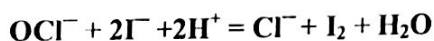


The liberated iodine is titrated with a standard solution of sodium thiosulphate.

**Procedure:**

1. Weigh out accurately about **2.5 gm** (weight taken = X) of the bleaching powder into a clean ceramic mortar. Add a little water, and rub the mixture to a smooth paste. Add a little more water, triturate with the **paste**, allow the mixture to settle and pour off the milky liquid into a **250 ml** volumetric flask. Grind the residue with a little more water and repeat the operation until the whole of the sample has been transferred to the flask either in solution or in a state of very fine **suspension**, and the mortar washed quite clean.
2. **Add water** to the volumetric flask up to mark and shake well.
3. Transfer **25 ml** of the turbid liquid immediately with a volumetric pipette into a **250ml** conical flask.
4. Add about **25 ml** of water into the conical flask.
5. Add about **2 gm** of solid **Iodate-free potassium iodide** to the conical flask.
6. Add **10 ml** of glacial acetic acid and shake, **Iodine** will be liberated immediately.
7. Keep the flask closed with a watch glass and place in a dark box for **5 minutes**. Then add **25 ml** distilled water to the flask.

8. Titrate the liberated iodine with standard sodium thiosulphate solution (standardized as per expt.04) until the color fades to pale yellow.
9. Add 1ml of starch indicator, a deep violet color appears.
10. Continue the titration with sodium thiosulphate until the color just disappears.

**Reactions:**

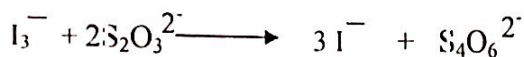

**Note:** Sodium thiosulphate should be standardized according to procedure (04) before/ after the experiment.

**Calculation:**

Available Chlorine content (%) in bleaching powder =

$$\frac{0.03546g \times a \times b \times V \times 100}{P \times Q}$$

Where **a** = volume of  $Na_2S_2O_3$  required in titration.

**b** = Molarity of  $Na_2S_2O_3$  solution.

**v** = Total volume of suspension (250ml).

**P** = Amount in gm of bleaching powder taken.

**Q** = Volume of suspension (25ml).taken for titration.

**Another Calculation:**

$$\text{Amount of bleaching powder in 50ml} = \frac{V \times E.O}{500} \text{ gm}$$

$$1\text{ml } 1M - Na_2S_2O_3 = 0.03546 \text{ g of Cl}$$

Experiment No. 10

**Name of the Experiment:** Drawing the pH-neutralization curves from the titration of a strong acid with a strong base and calculation the concentration of the strong acid.

Theory:

When alkali solution is gradually added to an acid solution, the pH of the solution increases due to neutralization of  $H^+$  ions. At the end point there is a sharp increase in the value of pH. The curve showing the variation of pH with the amount of alkali added is called the neutralization curve. The end point is marked by point of inflection in the curve. Further, knowledge of this curve allows selection of appropriate acid-base indicators for titration. For this purpose a pH meter will be used.

Procedure:

1. Take 10.0 ml of the supplied HCl solution in a 250 ml beaker, dilute to about 100 ml (90 ml water+10 ml HCl).
2. Place the electrode of the pH meter in the above solution.
3. Fill the burette with the supplied 0.1 M NaOH solution.
4. Add about 0.5 ml aliquot from the burette to the solution in the beaker. Stir the solution and measure its pH. After adding 9 ml NaOH solution then gradually add NaOH solution with 0.2 ml interval upto 11 ml. Add NaOH solution 1ml to 15 ml.
5. Tabulate burette readings, volume of alkali added, and pH of the solution.
6. Plot pH vs volume of NaOH solution added.
7. Draw a smooth curve; locate the equivalent point (end point). This is the volume of alkali required to neutralize the acid.
8. Calculate the concentration of HCl solution in terms of Molarity and Normality.

Questions:

1. What is meant by pH? Calculate the pH values of 0.1 M  $H_2SO_4$  and 0.1 M NaOH.
2. Draw pH – Neutralization curves for the titration of different types of acids and bases ( ex. strong acid and strong base) and comment on the suitable indicators for those titrations.



Calibration of pH meter

1. Wash the probe/ electrode of pH meter with distilled water.
2. Immerse the probe and temperature probe approximately 1.5 inch into a buffer solution having pH 7.0
3. Stir the beaker gently.
4. Press 'CAL'.
5. Not ready indication will blink on LCD until the ready was stabilized.
6. When the ready is stable and close to the selected buffer (indicated in the right-top side of the monitor) 'ready' indication will appear on 'CFM' will blink then press CFM to confirm the calibration. (it may take one or two minutes to stabilize the pH meter).
7. The calibrated value is then displayed on the right-top side of the LCD will display the second-expected buffer value.
8. Repeat the same procedure for the calibration of pH meter for pH 4.0.
9. Press CAL for quit calibration.


Experiment No. 11

**Name of the Experiment:** Determination of the strength of Hydrogen peroxide (H<sub>2</sub>O<sub>2</sub>)

Theory:

Hydrogen peroxide (H<sub>2</sub>O<sub>2</sub>), generally known as an oxidizer, is commonly used as a bleaching agent. It is the simplest peroxide (a compound with an oxygen-oxygen single bond). Hydrogen peroxide is a clear liquid, more viscous than water that appears colorless in dilute solution. It is also used as a disinfectant, antiseptic etc. The oxidizing capacity of hydrogen peroxide is so strong that it is considered a highly oxidizing species. Hydrogen peroxide is naturally produced in organisms as a by-product of oxidative metabolism.

Though H<sub>2</sub>O<sub>2</sub> usually act as an oxidizing agent (oxidizer), in the presence of a stronger oxidizing agent like KMnO<sub>4</sub>, it behaves as a reducing agent. When potassium permanganate solution (KMnO<sub>4</sub>) is added to hydrogen peroxide (H<sub>2</sub>O<sub>2</sub>) solution in the presence of sulphuric acid the following reaction occurs. This principle is used to determine the strength of the H<sub>2</sub>O<sub>2</sub>.

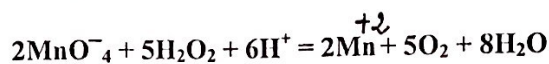

Procedure:

1. Transfer 10 ml of the supplied solution of H<sub>2</sub>O<sub>2</sub> into a conical flask.
2. Dilute this solution with 25 ml of water and shake thoroughly.
3. Add 10 ml of 1M H<sub>2</sub>SO<sub>4</sub>.
4. Perform the titration with standard 0.1M potassium permanganate to the first permanent, faint pink color.
5. Repeat the titration three times.

**Note:** Standardize the potassium permanganate solution (exp: 06).

Calculation:

**Method 1**

Find out the concentration of H<sub>2</sub>O<sub>2</sub>. Then calculate the % as instructed by teacher.

**Method 2**

1. The weight of hydrogen peroxide per litre of the original solution.
2. The 'volume strength' i.e the number of ml of oxygen at N.T.P that can be obtained from 1ml of the original solution.

$$1 \text{ ml } 1 \text{ M KMnO}_4 \text{ solution} = 0.01701 \times 5 \text{ g H}_2\text{O}_2$$

$$= 0.08505 \text{ g H}_2\text{O}_2$$

$$\text{Concentration of H}_2\text{O}_2 = \frac{0.08505 \times \text{vol. of KMnO}_4 \text{ in ml} \times \text{Molarity of KMnO}_4}{10} \times 100$$

(% by volume)

**Questions:**

1. Give half reactions and explain how  $\text{H}_2\text{O}_2$  act as a reducing agent.
2. What is bleaching?




Experiment No. 12

Name of the Experiment: Preparation of Iron (II) Sulphate heptahydrate  $\text{FeSO}_4 \cdot 7\text{H}_2\text{O}$

Theory:

When metallic iron (Fe) reacts with sulphuric acid, it replaces the acidic hydrogen ( $\text{H}^+$ ) in  $\text{H}_2\text{SO}_4$  and produces iron (II) sulphate.

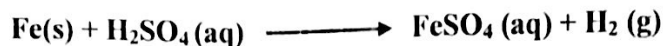

$\text{FeSO}_4$  crystallizes from the aqueous solution with seven water molecules.

Procedure:

1. Pour **10 ml** of concentrated sulphuric acid slowly into **200 ml** of distilled water in a 500 ml beaker and add **10 gm** of iron turnings.
2. When the reaction **slows** down, heat the beaker until the acid is practically all neutralized as indicated by the fact that evolution of hydrogen gas ceases.
3. **Filter** the content of the beaker, using a filter paper to remove any undissolved residue.
4. The filtrate should give a **bottle green** color.
5. If the solution is oxidized, appreciable amount of sulphuric acid must be added until the color becomes green. Oxidation will be indicated by the change in color of the solution from bottle green to **yellowish shade** of green or by the formation of rusty precipitate.
6. **Slowly** evaporate the solution to crystallization and cool.
7. Then place the beaker in an ice bath to complete crystallization.
8. **Filter off** the crystals from the mother liquor.
9. **Wash** the crystals twice with small portions of mother liquor.
10. **Dry** the crystals using absorbent paper.
11. **Weigh** the product and determine the yield.
12. Store the crystal in sample bottle.

Calculation:

1. Calculate the no. of moles of iron you have taken in this experiment.
2. Calculate the no. of moles of iron (II) sulphate that you have obtained in this experiment.

Questions:


3. What happens when iron (II) sulphate is kept in contact with air?
4. What type of reaction takes place during the conversion of iron to ferrous sulphate?
5. Write down the electronic configuration of iron (II) and iron (III).

Experiment No.s13

Name of the Experiment: Determination of  $\text{Na}_2\text{CO}_3$  content in commercial Soda.

Theory:

Sodium carbonate ( $\text{Na}_2\text{CO}_3$ ) content in commercial caustic soda solution can be determined by titrimetric method. The alkali carbonate ( $\text{CO}_3^{2-}$ ) of the supplied solution can be determined by titrating with a standard acid solution. Then the carbonate ion in the solution may be precipitated by a dilute solution of  $\text{BaCl}_2$ .

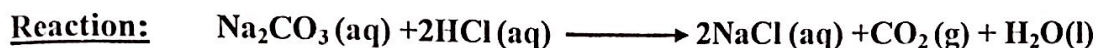

Procedure:

1. Weigh out accurately 1.44 gm of washing soda in a 100 ml volumetric flask.
2. Dissolve with a little water and make up to mark.
3. Take 10 ml of the solution with a pipette in a 250 ml conical flask.
4. Add about 25 ml of water.
5. Add 1-2 drops of methyl orange indicator to the solution.
6. Titrate the solution with standard hydrochloric acid solution (standardized as per expr:03) until the color changes from yellow to orange.
7. Repeat the procedure at three at least three times and note the mean value.

Calculation:

1 ml N HCl = 0.05300 gm  $\text{Na}_2\text{CO}_3$

$$\% \text{ of } \text{Na}_2\text{CO}_3 \text{ in washing soda} = \frac{0.05300 \text{ gm} \times a \times b \times c}{P \times Q} \times 100$$

Where a = volume of HCl (ml) solution.

b = Molarity of HCl solution.

c = Total volume of solution (100ml).

P = Volume taken for titration (10ml).

Q = Weight of soda in gm.

### Report writing

- ♣ Report on each experiment should be submitted on the next class.
- ♣ Report should be submitted separately for each experiment.
- ♣ A4 size offset paper should be used. Writing must be on one side.

#### ♣ **Top sheet of the report must contain the following points**

Chemistry sessional report, Course no., Experiment no, Name of the experiment, Date of Performance, Date of Submission, Performer's Name, Department, ID, Group no., Group Partner's ID.

- ♣ Top sheet may be composed but the rest of the report must be hand written.

#### ▶ Report should contain the followings:

- ✍ Name of the experiment
- ✍ Objective of the experiment
- Theory (related to the experiment, equation, half-reactions etc.)
- ✍ Apparatus required
- ✍ Chemicals required
- ✍ Preparation of standard solution (passive form to be used) with calculation.
- ✍ Experimental data table with title.
- ✍ Calculation (both in formula and unitary methods)
- Result
- Major cautions if any
- Answers to the questions.

### Students Must Learn

Qualitative analysis and Quantitative analysis

#### *Titration*

Primary standard substance, Secondary standard substance

Indicator, Auto indicator, End point

*Solution*, Solute, Solvent, Dilute solution, Standard solution

#### *Concentration*

Expression of Concentration: Molarity, Molality, Normality, %, mole fraction



*Neutralization reaction*

Acids and Bases (**Arrhenius**, Brönsted & Lewis concepts), Salt  
pH, Buffer solution, Titration curves - suitable indicator, Theories of indicator

*Redox reaction*

Oxidation, Reduction half-reactions, Oxidizing agent, Reducing agent

Iodometry and Iodimetry

*Hardness of water*

Ligand, Complex, Chelating agent, Chelate

**DOs and DON'Ts (Cautions)**

- All students must wear apron (long, full sleeve, white) during the lab work.
- Students must not do anything without any given instruction.
- Students must handle apparatus and chemicals carefully.
- Put on Safety glass.
- Maintain cleanliness.
- Be disciplined.

**References**

1. Essentials of Physical Chemistry – by Bahl & Tuli
2. General Chemistry – by Ebbing
3. Modern Inorganic Chemistry – by S. Z. Haider
4. Vogel's Text Book of Quantitative Chemical Analysis (Pearson Education Asia)
4. Analytical Chemistry – by Gary D. Christian (John Wiley & Sons Inc.)